


Advazorb

Hydrophilic Foam Dressing Range

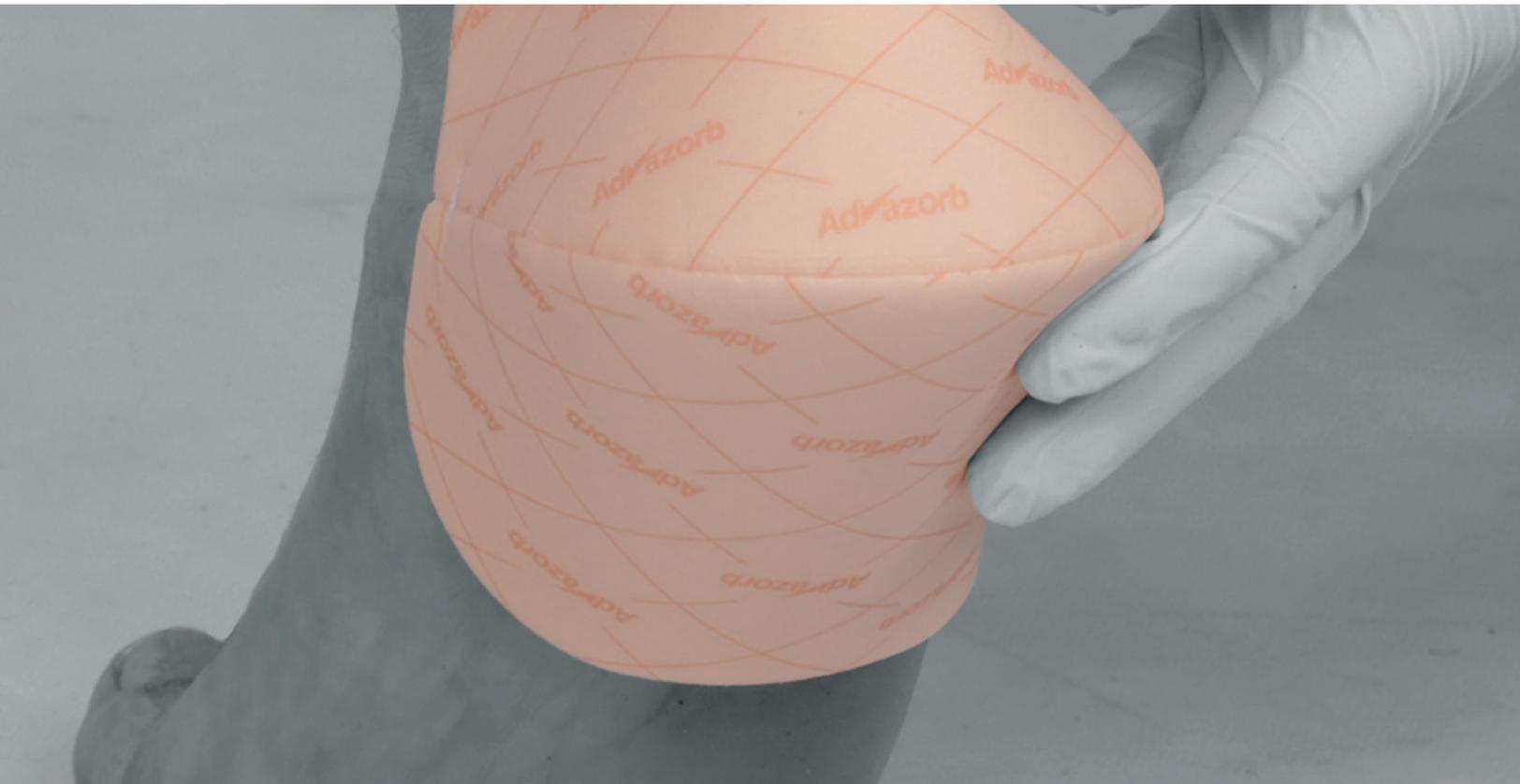

Gentle on skin

Distributed by



PASSION. PARTNERSHIP. POSSIBILITIES.

www.dukal.com



www.advancismedical.com


Advazorb®

Hydrophilic foam dressing range

Manufactured in the UK, the Advancis Medical® range of Advazorb® Hydrophilic foam dressings, offer the combination of a high quality performance product at a competitive price.

Advazorb® is a comprehensive range of patient friendly, absorbent foam dressings presented in non-adhesive and atraumatic silicone adhesive format. The range is designed to overcome the complex challenges of managing exudate whilst protecting 'at risk' fragile skin, a combination of excellent fluid handling capabilities of Advazorb® with a layer of atraumatic soft silicone adhesive.

Advazorb®, Advazorb Silfix® and Advazorb Border® are available in both a regular thickness for moderate to high exudate and 'Life' versions for low to moderately exuding wounds.

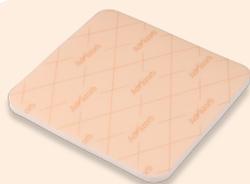

Advazorb®
Hydrophilic foam dressing
with film backing

Soft and conformable, non-adherent, hydrophilic, polyurethane foam dressings with a breathable film backing.

A non-adherent foam dressing with excellent fluid handling properties. Ideal for use under compression bandaging with its fluid retention, low profile and low friction film backing.

- High quality polyurethane foam, rapidly absorbs and retains what it absorbs. Available in regular and 'Life' versions to suit a wider range of wounds.
- The Advazorb® range features a waterproof film backing that prevents strike-through. The backing has been designed to reduce friction on clothes and bandages which extends dressing wear time by reducing 'rucking' and preventing dressings falling off.
- With an MVTR range of 2600-2700 g/m²/24hours the Advazorb® foam range provides the ideal conditions for wound healing.
- Increasing patient compliance and comfort.
- Can be used on a wide variety of exuding wounds including under compression and can be cut to shape and size.



Advazorb Heel®
Hydrophilic foam dressing
with film backing

A conformable, non-adherent, hydrophilic, polyurethane foam dressing with a breathable film backing designed to fit around the heel.

A non-adherent foam dressing with excellent fluid handling properties. Ideal for use under compression bandaging with its fluid retention, low profile and low friction film backing.

- The same high quality PU foam as the other dressings in the Advazorb® range, it has excellent fluid handling properties and rapidly absorbs and retains exudate.
- It features the low friction, bacteria-proof, thin film backing that the entire Advazorb® range has meaning extended wear time, reduced 'rucking' and is waterproof.
- Soft and conformable, increasing patient comfort and is ideal to use under compression due to its low profile.
- It is anatomically shaped to fit around the heel so that the dressing is in close contact with the wound to provide an environment conducive to healing.

Item #	Description	Packaging	HCPCS
CR4165	Advazorb®, 3" x 3"	10/bx	A6209
CR4166	Advazorb®, 4" x 4"	10/bx	A6209
CR4167	Advazorb®, 5" x 5"	10/bx	A6210
CR4168	Advazorb®, 6" x 6"	10/bx	A6210
CR4169	Advazorb®, 4" x 8"	10/bx	A6210
CR4170	Advazorb®, 8" x 8"	10/bx	A6211
CR4228	Advazorb Heel®, 6.5" x 8.5"	5/bx	A6210